



# Cambridge International AS & A Level

CANDIDATE  
NAME

CENTRE  
NUMBER

--	--	--	--	--	--

CANDIDATE  
NUMBER

--	--	--	--

\* 5 6 1 6 0 2 1 7 0 8 \*



**CHEMISTRY**

**9701/35**

Paper 3 Advanced Practical Skills 1

**May/June 2021**

**2 hours**

You must answer on the question paper.

You will need: The materials and apparatus listed in the confidential instructions

## INSTRUCTIONS

- Answer **all** questions.
- Use a black or dark blue pen. You may use an HB pencil for any diagrams or graphs.
- Write your name, centre number and candidate number in the boxes at the top of the page.
- Write your answer to each question in the space provided.
- Do **not** use an erasable pen or correction fluid.
- Do **not** write on any bar codes.
- You may use a calculator.
- You should show all your working, use appropriate units and use an appropriate number of significant figures.
- Give details of the practical session and laboratory, where appropriate, in the boxes provided.

<b>Session</b>	
<b>Laboratory</b>	

## INFORMATION

- The total mark for this paper is 40.
- The number of marks for each question or part question is shown in brackets [ ].
- The Periodic Table is printed in the question paper.
- Notes for use in qualitative analysis are provided in the question paper.

<b>For Examiner's Use</b>	
<b>1</b>	
<b>2</b>	
<b>3</b>	
<b>Total</b>	

This document has **12** pages. Any blank pages are indicated.

## Quantitative analysis

Read through the whole method before starting any practical work. Where appropriate, prepare a table for your results in the space provided.

Show your working and appropriate significant figures in the final answer to **each** step of your calculations.

- 1 In this experiment you will carry out a titration to determine the relative formula mass,  $M_r$ , of a monoprotic acid, HX.

**FA 1** is HX, the monoprotic acid.

**FA 2** is  $0.0450 \text{ mol dm}^{-3}$  sodium carbonate,  $\text{Na}_2\text{CO}_3$ .  
methyl orange indicator

### (a) Method

#### Preparing a solution of FA 1

- Weigh the empty  $250 \text{ cm}^3$  beaker. Record the mass.
- Transfer all the **FA 1** into the beaker.
- Weigh the beaker and **FA 1**. Record the mass.
- Calculate and record the mass of **FA 1** used.
- Add approximately  $100 \text{ cm}^3$  of distilled water to the **FA 1** in the beaker.
- Stir the mixture with a glass rod until all the **FA 1** has dissolved.
- Transfer this solution into the  $250 \text{ cm}^3$  volumetric flask.
- Wash the beaker with distilled water and transfer the washings to the volumetric flask.
- Rinse the glass rod with distilled water and transfer the washings to the volumetric flask.
- Make up the solution in the volumetric flask to the mark using distilled water.
- Shake the flask thoroughly.
- This solution of HX is **FA 3**. Label the flask **FA 3**.

### Results

#### Titration

- Fill the burette with **FA 3**.
- Pipette  $25.0 \text{ cm}^3$  of **FA 2** into a conical flask.
- Add several drops of methyl orange indicator.
- Perform a rough titration and record your burette readings in the space below.

The rough titre is .....  $\text{cm}^3$ .

- Carry out as many accurate titrations as you think necessary to obtain consistent results.
- Make sure any recorded results show the precision of your practical work.
- Record in a suitable form below all of your burette readings and the volume of **FA 3** added in each accurate titration.

I	
II	
III	
IV	
V	
VI	
VII	
VIII	

[8]

- (b) From your accurate titration results, obtain a suitable value for the volume of **FA 3** to be used in your calculations.  
Show clearly how you obtained this value.


25.0 cm<sup>3</sup> of **FA 2** required ..... cm<sup>3</sup> of **FA 3**. [1]

**(c) Calculations**

- (i) Calculate the number of moles of sodium carbonate present in the volume of **FA 2** used in each titration.

moles of Na<sub>2</sub>CO<sub>3</sub> = ..... mol [1]

- (ii) Give the equation for the reaction of **FA 2**, Na<sub>2</sub>CO<sub>3</sub>, with **FA 3**, HX.

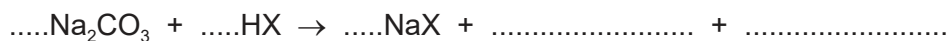

Use your answer to **(c)(i)** to deduce the number of moles of HX present in the volume you calculated in **(b)**.

moles of HX = ..... mol [1]

- (iii) Use your answer to (c)(ii) and your data recorded on page 2 to calculate the relative formula mass,  $M_r$ , of HX.  
Show your working.

$M_r$  of HX = ..... [2]

- (iv) One molecule of HX contains one nitrogen atom, three oxygen atoms, three hydrogen atoms and one atom of another element, **E**.  
The identity of **E** can be found by calculation.

Show this calculation and identify **E**.

Element **E** is ..... [1]

- (d) What is the error in a single reading for the balance that you used?

error = ..... g

Calculate the maximum percentage error in the mass of **FA 1** that you recorded on page 2.

maximum percentage error = ..... %  
[1]

- (e) Suggest and carry out an experiment using aqueous silver nitrate to determine whether the compound AgX is soluble or insoluble in water.

method .....

.....

observations .....

conclusion .....

[2]

- (f) Suggest why the use of methyl orange indicator might give an inaccurate titration result.

.....

..... [1]

[Total: 18]

- 2 You will determine the percentage by mass of water of crystallisation in a hydrated salt by thermal decomposition. When a hydrated salt is heated it loses its water of crystallisation.

**FA 4** is the hydrated salt.

**(a) Method**

- Weigh the crucible with its lid. Record the mass.
- Transfer 2.0–2.2 g of **FA 4** from the container into the crucible.
- Weigh the crucible, lid and **FA 4**. Record the mass.
- Calculate and record the mass of **FA 4** used.
- Place the crucible and contents on the pipe-clay triangle.
- Heat gently, with the lid on, for approximately 1 minute.
- Heat strongly, with the lid off, for a further 4 minutes.
- Replace the lid and leave the crucible to cool for at least 5 minutes.

**During the cooling period, you may wish to begin work on Question 3.**

- When the crucible is cool, weigh the crucible with its lid and contents. Record the mass.
- Heat strongly, with the lid off, for a further 2 minutes.
- Replace the lid and leave the crucible to cool for at least 5 minutes.
- When it is cool, weigh the crucible with its lid and contents. Record the mass.
- Calculate and record the mass lost.

**Keep FA 4 for use in Question 3.**

**Results**

I	
II	
III	
IV	
V	

[5]

- (b) (i)** Calculate the percentage by mass of water of crystallisation in **FA 4**.

percentage by mass of H<sub>2</sub>O in **FA 4** = ..... % [1]

- (ii)** State what assumption you made about anhydrous **FA 4** in your calculation in **(b)(i)**.

.....  
 ..... [1]

[Total: 7]

## Qualitative analysis

Where reagents are selected for use in a test, the **name** or **correct formula** of the element or compound must be given.

At each stage of any test you are to record details of the following:

- colour changes seen
- the formation of any precipitate and its solubility in an excess of the reagent added
- the formation of any gas and its identification by a suitable test.

You should indicate clearly at what stage in a test a change occurs.

If any solution is warmed, a **boiling tube** must be used.

Rinse and reuse test-tubes and boiling tubes where possible.

**No additional tests for ions present should be attempted.**

- 3 (a) **FA 4** contains one cation and one anion both of which are listed in the Qualitative Analysis Notes.

Dissolve the remainder of the **FA 4** in approximately 20 cm<sup>3</sup> of water in the 100 cm<sup>3</sup> beaker. Use a 1 cm depth of this solution in a test-tube for each of the tests.

- (i) The anion in **FA 4** is chloride.

State which reagents you would use to identify this ion conclusively.

.....

Use your reagents to test your solution of **FA 4**.  
Record your observations.

.....

.....

.....

In what way are your observations **not** consistent with the expected results for chloride?

.....

.....

[3]

- (ii) Add an equal volume of dilute sulfuric acid to aqueous **FA 4**.

Record your observations and deduce the identity of the cation in **FA 4**.

observations .....

identity of cation .....

Give the ionic equation for this reaction. Include state symbols.

..... [2]

- (b) Use the conclusions in (a)(i) and (a)(ii) to calculate the relative formula mass of anhydrous **FA 4**.

$M_r =$  .....

Use your results from **Question 2** to calculate the number of moles of water of crystallisation in one mole of hydrated **FA 4**.

moles of  $H_2O =$  ..... [3]

- (c) **FA 5** contains one cation and one anion, both of which are listed in the Qualitative Analysis Notes.

- (i) Heat a small spatula measure of **FA 5** in a hard-glass test-tube. Keep heating until no further change is observed. Record **all** your observations.

.....  
 .....  
 .....  
 ..... [2]

- (ii) To a 1 cm depth of dilute sulfuric acid in a test-tube, add a small spatula measure of **FA 5**. Record all your observations.

.....

.....

.....

.....

..... [3]

- (iii) Place a small spatula measure of **FA 5** into a boiling tube. Add a 1 cm depth of aqueous sodium hydroxide and warm the mixture gently. Record your observations.

.....

..... [1]

- (iv) From the results of your tests, deduce the chemical formula of **FA 5**.

**FA 5** is ..... [1]

[Total: 15]





## Qualitative analysis notes

## 1 Reactions of aqueous cations

ion	reaction with	
	NaOH(aq)	NH <sub>3</sub> (aq)
aluminium, Al <sup>3+</sup> (aq)	white ppt. soluble in excess	white ppt. insoluble in excess
ammonium, NH <sub>4</sub> <sup>+</sup> (aq)	no ppt. ammonia produced on heating	–
barium, Ba <sup>2+</sup> (aq)	faint white ppt. is nearly always observed unless reagents are pure	no ppt.
calcium, Ca <sup>2+</sup> (aq)	white ppt. with high [Ca <sup>2+</sup> (aq)]	no ppt.
chromium(III), Cr <sup>3+</sup> (aq)	grey-green ppt. soluble in excess	grey-green ppt. insoluble in excess
copper(II), Cu <sup>2+</sup> (aq)	pale blue ppt. insoluble in excess	pale blue ppt. soluble in excess giving dark blue solution
iron(II), Fe <sup>2+</sup> (aq)	green ppt. turning brown on contact with air insoluble in excess	green ppt. turning brown on contact with air insoluble in excess
iron(III), Fe <sup>3+</sup> (aq)	red-brown ppt. insoluble in excess	red-brown ppt. insoluble in excess
magnesium, Mg <sup>2+</sup> (aq)	white ppt. insoluble in excess	white ppt. insoluble in excess
manganese(II), Mn <sup>2+</sup> (aq)	off-white ppt. rapidly turning brown on contact with air insoluble in excess	off-white ppt. rapidly turning brown on contact with air insoluble in excess
zinc, Zn <sup>2+</sup> (aq)	white ppt. soluble in excess	white ppt. soluble in excess

## 2 Reactions of anions

<i>ion</i>	<i>reaction</i>
carbonate, $\text{CO}_3^{2-}$	$\text{CO}_2$ liberated by dilute acids
chloride, $\text{Cl}^-(\text{aq})$	gives white ppt. with $\text{Ag}^+(\text{aq})$ (soluble in $\text{NH}_3(\text{aq})$ )
bromide, $\text{Br}^-(\text{aq})$	gives cream ppt. with $\text{Ag}^+(\text{aq})$ (partially soluble in $\text{NH}_3(\text{aq})$ )
iodide, $\text{I}^-(\text{aq})$	gives yellow ppt. with $\text{Ag}^+(\text{aq})$ (insoluble in $\text{NH}_3(\text{aq})$ )
nitrate, $\text{NO}_3^-(\text{aq})$	$\text{NH}_3$ liberated on heating with $\text{OH}^-(\text{aq})$ and $\text{Al}$ foil
nitrite, $\text{NO}_2^-(\text{aq})$	$\text{NH}_3$ liberated on heating with $\text{OH}^-(\text{aq})$ and $\text{Al}$ foil
sulfate, $\text{SO}_4^{2-}(\text{aq})$	gives white ppt. with $\text{Ba}^{2+}(\text{aq})$ (insoluble in excess dilute strong acids)
sulfite, $\text{SO}_3^{2-}(\text{aq})$	gives white ppt. with $\text{Ba}^{2+}(\text{aq})$ (soluble in excess dilute strong acids)

## 3 Tests for gases

<i>gas</i>	<i>test and test result</i>
ammonia, $\text{NH}_3$	turns damp red litmus paper blue
carbon dioxide, $\text{CO}_2$	gives a white ppt. with limewater (ppt. dissolves with excess $\text{CO}_2$ )
chlorine, $\text{Cl}_2$	bleaches damp litmus paper
hydrogen, $\text{H}_2$	'pops' with a lighted splint
oxygen, $\text{O}_2$	relights a glowing splint

## The Periodic Table of Elements

		Group															
1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18
		<div style="display: flex; justify-content: space-around; align-items: center;"> <div style="border: 1px solid black; padding: 2px;">1 H hydrogen 1.0</div> <div style="border: 1px solid black; padding: 2px;">2 He helium 4.0</div> </div>															
		<div style="display: flex; justify-content: space-around; align-items: center;"> <div style="border: 1px solid black; padding: 2px;">3 Li lithium 6.9</div> <div style="border: 1px solid black; padding: 2px;">4 Be beryllium 9.0</div> <div style="border: 1px solid black; padding: 2px;">5 B boron 10.8</div> <div style="border: 1px solid black; padding: 2px;">6 C carbon 12.0</div> <div style="border: 1px solid black; padding: 2px;">7 N nitrogen 14.0</div> <div style="border: 1px solid black; padding: 2px;">8 O oxygen 16.0</div> <div style="border: 1px solid black; padding: 2px;">9 F fluorine 19.0</div> <div style="border: 1px solid black; padding: 2px;">10 Ne neon 20.2</div> </div>															
		<div style="display: flex; justify-content: space-around; align-items: center;"> <div style="border: 1px solid black; padding: 2px;">11 Na sodium 23.0</div> <div style="border: 1px solid black; padding: 2px;">12 Mg magnesium 24.3</div> <div style="border: 1px solid black; padding: 2px;">13 Al aluminium 27.0</div> <div style="border: 1px solid black; padding: 2px;">14 Si silicon 28.1</div> <div style="border: 1px solid black; padding: 2px;">15 P phosphorus 31.0</div> <div style="border: 1px solid black; padding: 2px;">16 S sulfur 32.1</div> <div style="border: 1px solid black; padding: 2px;">17 Cl chlorine 35.5</div> <div style="border: 1px solid black; padding: 2px;">18 Ar argon 39.9</div> </div>															
		<div style="display: flex; justify-content: space-around; align-items: center;"> <div style="border: 1px solid black; padding: 2px;">19 K potassium 39.1</div> <div style="border: 1px solid black; padding: 2px;">20 Ca calcium 40.1</div> <div style="border: 1px solid black; padding: 2px;">21 Sc scandium 45.0</div> <div style="border: 1px solid black; padding: 2px;">22 Ti titanium 47.9</div> <div style="border: 1px solid black; padding: 2px;">23 V vanadium 50.9</div> <div style="border: 1px solid black; padding: 2px;">24 Cr chromium 52.0</div> <div style="border: 1px solid black; padding: 2px;">25 Mn manganese 54.9</div> <div style="border: 1px solid black; padding: 2px;">26 Fe iron 55.8</div> <div style="border: 1px solid black; padding: 2px;">27 Co cobalt 58.9</div> <div style="border: 1px solid black; padding: 2px;">28 Ni nickel 58.7</div> <div style="border: 1px solid black; padding: 2px;">29 Cu copper 63.5</div> <div style="border: 1px solid black; padding: 2px;">30 Zn zinc 65.4</div> <div style="border: 1px solid black; padding: 2px;">31 Ga gallium 69.7</div> <div style="border: 1px solid black; padding: 2px;">32 Ge germanium 72.6</div> <div style="border: 1px solid black; padding: 2px;">33 As arsenic 74.9</div> <div style="border: 1px solid black; padding: 2px;">34 Se selenium 79.0</div> <div style="border: 1px solid black; padding: 2px;">35 Br bromine 79.9</div> <div style="border: 1px solid black; padding: 2px;">36 Kr krypton 83.8</div> </div>															
		<div style="display: flex; justify-content: space-around; align-items: center;"> <div style="border: 1px solid black; padding: 2px;">37 Rb rubidium 85.5</div> <div style="border: 1px solid black; padding: 2px;">38 Sr strontium 87.6</div> <div style="border: 1px solid black; padding: 2px;">39 Y yttrium 88.9</div> <div style="border: 1px solid black; padding: 2px;">40 Zr zirconium 91.2</div> <div style="border: 1px solid black; padding: 2px;">41 Nb niobium 92.9</div> <div style="border: 1px solid black; padding: 2px;">42 Mo molybdenum 95.9</div> <div style="border: 1px solid black; padding: 2px;">43 Tc technetium —</div> <div style="border: 1px solid black; padding: 2px;">44 Ru ruthenium 101.1</div> <div style="border: 1px solid black; padding: 2px;">45 Rh rhodium 102.9</div> <div style="border: 1px solid black; padding: 2px;">46 Pd palladium 106.4</div> <div style="border: 1px solid black; padding: 2px;">47 Ag silver 107.9</div> <div style="border: 1px solid black; padding: 2px;">48 Cd cadmium 112.4</div> <div style="border: 1px solid black; padding: 2px;">49 In indium 114.8</div> <div style="border: 1px solid black; padding: 2px;">50 Sn tin 118.7</div> <div style="border: 1px solid black; padding: 2px;">51 Sb antimony 121.8</div> <div style="border: 1px solid black; padding: 2px;">52 Te tellurium 127.6</div> <div style="border: 1px solid black; padding: 2px;">53 I iodine 126.9</div> <div style="border: 1px solid black; padding: 2px;">54 Xe xenon 131.3</div> </div>															
		<div style="display: flex; justify-content: space-around; align-items: center;"> <div style="border: 1px solid black; padding: 2px;">55 Cs caesium 132.9</div> <div style="border: 1px solid black; padding: 2px;">56 Ba barium 137.3</div> <div style="border: 1px solid black; padding: 2px;">57–71 lanthanoids</div> <div style="border: 1px solid black; padding: 2px;">72 Hf hafnium 178.5</div> <div style="border: 1px solid black; padding: 2px;">73 Ta tantalum 180.9</div> <div style="border: 1px solid black; padding: 2px;">74 W tungsten 183.8</div> <div style="border: 1px solid black; padding: 2px;">75 Re rhenium 186.2</div> <div style="border: 1px solid black; padding: 2px;">76 Os osmium 190.2</div> <div style="border: 1px solid black; padding: 2px;">77 Ir iridium 192.2</div> <div style="border: 1px solid black; padding: 2px;">78 Pt platinum 195.1</div> <div style="border: 1px solid black; padding: 2px;">79 Au gold 197.0</div> <div style="border: 1px solid black; padding: 2px;">80 Hg mercury 200.6</div> <div style="border: 1px solid black; padding: 2px;">81 Tl thallium 204.4</div> <div style="border: 1px solid black; padding: 2px;">82 Pb lead 207.2</div> <div style="border: 1px solid black; padding: 2px;">83 Bi bismuth 209.0</div> <div style="border: 1px solid black; padding: 2px;">84 Po polonium —</div> <div style="border: 1px solid black; padding: 2px;">85 At astatine —</div> <div style="border: 1px solid black; padding: 2px;">86 Rn radon —</div> </div>															
		<div style="display: flex; justify-content: space-around; align-items: center;"> <div style="border: 1px solid black; padding: 2px;">87 Fr francium —</div> <div style="border: 1px solid black; padding: 2px;">88 Ra radium —</div> <div style="border: 1px solid black; padding: 2px;">89–103 actinoids</div> <div style="border: 1px solid black; padding: 2px;">104 Rf rutherfordium —</div> <div style="border: 1px solid black; padding: 2px;">105 Db dubnium —</div> <div style="border: 1px solid black; padding: 2px;">106 Sg seaborgium —</div> <div style="border: 1px solid black; padding: 2px;">107 Bh bohrium —</div> <div style="border: 1px solid black; padding: 2px;">108 Hs hassium —</div> <div style="border: 1px solid black; padding: 2px;">109 Mt meitnerium —</div> <div style="border: 1px solid black; padding: 2px;">110 Ds darmstadtium —</div> <div style="border: 1px solid black; padding: 2px;">111 Rg roentgenium —</div> <div style="border: 1px solid black; padding: 2px;">112 Cn copernicium —</div> <div style="border: 1px solid black; padding: 2px;">114 Fl flerovium —</div> <div style="border: 1px solid black; padding: 2px;">116 Lv livermorium —</div> </div>															

lanthanoids

57	La	lanthanum	138.9	58	Ce	cerium	140.1	59	Pr	praseodymium	140.9	60	Nd	neodymium	144.4	61	Pm	promethium	—	62	Sm	samarium	150.4	63	Eu	europlium	152.0	64	Gd	gadolinium	157.3	65	Tb	terbium	158.9	66	Dy	dysprosium	162.5	67	Ho	holmium	164.9	68	Er	erbium	167.3	69	Tm	thulium	168.9	70	Yb	ytterbium	173.1	71	Lu	lutetium	175.0
89	Ac	actinium	—	90	Th	thorium	232.0	91	Pa	protactinium	231.0	92	U	uranium	238.0	93	Np	neptunium	—	94	Pu	plutonium	—	95	Am	americium	—	96	Cm	curium	—	97	Bk	berkelium	—	98	Cf	californium	—	99	Es	einsteinium	—	100	Fm	fermium	—	101	Md	mendeleevium	—	102	No	nobelium	—	103	Lr	lawrencium	—

actinoids